



Multi-Directional headrest support



PHYSIPRO

Seating and mobility solutions



This headrest support made of aluminum offers lightness, safety and allows adjustment in a variety of positions. With continuous adjustment, it can be deported in depth.

Adjustments

9" horizontal adjustment from front to back

6½" vertical adjustment

45° adjustment to the left and right



Headrest options

Basic headrest pad



Combination head/neck support pad



Adjustable lateral control headrest pad - Single



Adjustable lateral control headrest pad

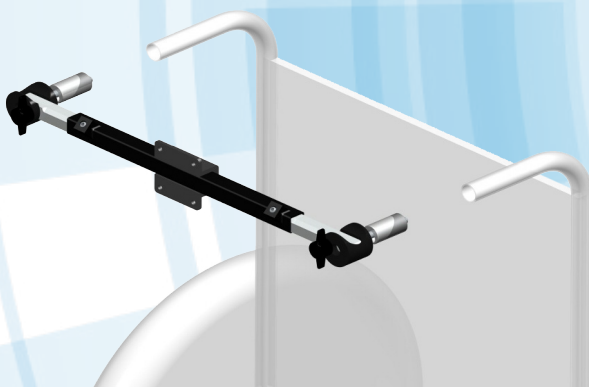


Headrest support mounting fixture

Internal

This swing-away and retractable system provides greater firmness for flexible backrests by maintaining the backrest canes in their original position. The fixture is installed in the tubes of the backrest handles.


Adjustment range from :
13" to 20".



External





Installation and adjustments are simple with this easily retractable system used for mounting a headrest support on a flexible backrest, allowing the backrest canes to maintain their original position.

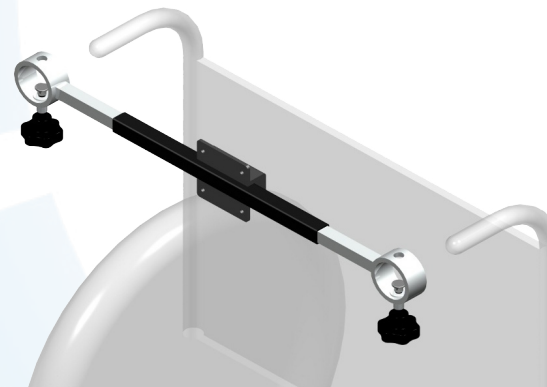
Adjustment range from :
13" to 18".





Technical specifications

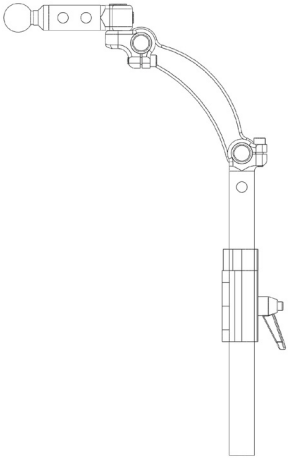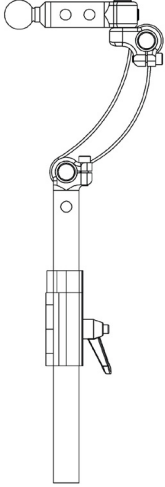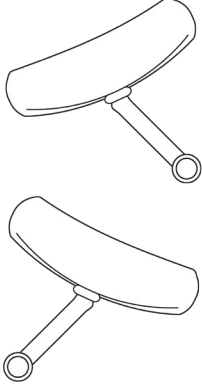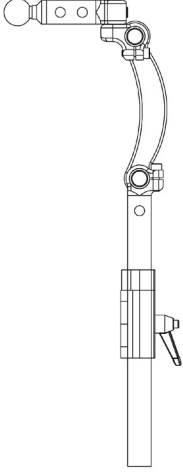Vertical adjustment	Horizontal adjustment	Horizontal axis	Total height
			
6½"	9"	45°	16"

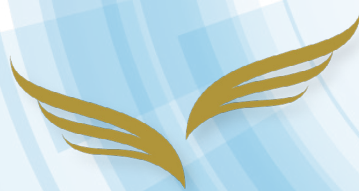

PHYSIPRO

Seating and mobility solutions

Head Office

370, 10th South Ave.,
Sherbrooke, Quebec J1G 2R7
T. 819 823-2252 | 1 800 668-2252
F. 819 565-3337
info@physipro.com
www.physipro.com

European branch

Village des Entrepreneurs
461, Saint-Léonard
49000 Angers, France
T. 02 41 69 38 01
F. 02 41 69 43 32
contact@physipro.fr