


StimuLITE

SILVER CUSHION



STIMULITE: SILVER CUSHION

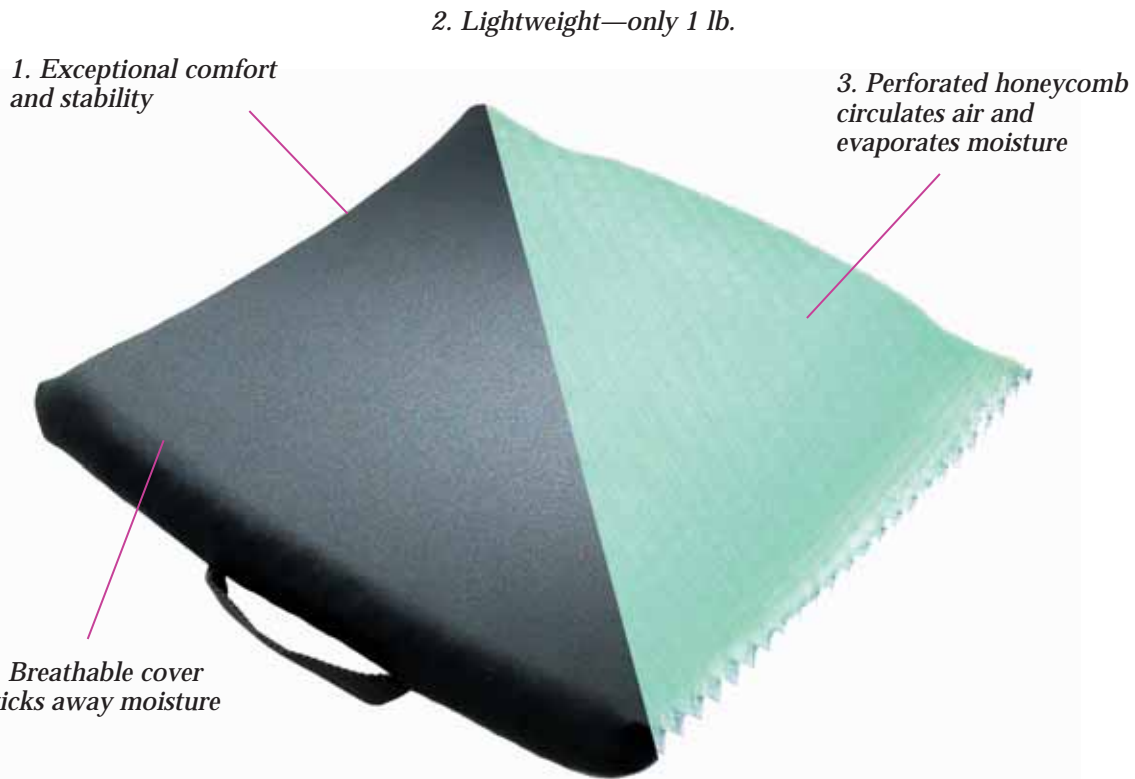

TOTAL PRESSURE[®] MANAGEMENT

Designed for elder care

This affordable, comfortable cushion can be used in a wheelchair, on a recliner, in a car, plane or on a stadium seat and provides exceptional stability and support. Made from a soft, flexible form of aerospace honeycomb, the lightweight Stimulite Silver Cushion distributes weight uniformly to relieve pressure. The honeycomb cells contour to body shape and flex with movement to absorb shock and prevent tissue trauma. Perforations in the honeycomb cell walls circulate air and evaporate moisture so you stay cool and comfortable, especially in hot and humid climates.



Made from an advanced medical grade material that's antibacterial, antifungal and odor resistant, this machine-washable cushion helps promote clean, healthy skin and is a breeze to care for.

Ideal for traveling, the Silver Cushion can be easily tossed into a car's trunk or an airplane's overhead compartment and washed in your hotel room's bathtub or shower. Just sponge with soapy water, rinse and shake to dry. A special breathable cover that wicks away moisture completes the honeycomb's unique ventilation system.




STIMULITE: SILVER CUSHION

Total Pressure Management



Maceration Resistance


Skin maceration is softening of the skin as a result of continual exposure to heat and bodily fluids or moisture. Severely softened or macerated skin is less tolerant of pressure and can lead to breakdown. Stimulite Honeycomb is perforated, enabling air to circulate and moisture to evaporate which helps maintain skin integrity. All of our Stimulite products are *machine washable* and *dryer safe*, making it easy to maintain good skin hygiene.

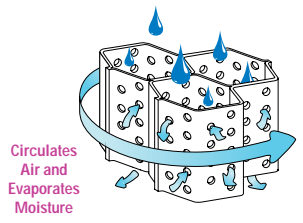

Complete Ventilation System

Pressure Relief plus Positioning


Honeycomb's cellular structure provides uniform load distribution and extraordinary comfort. Stimulite Honeycomb relieves pressure by distributing it away from hot spots. As a result, bony protuberances "feel" the same pressure as the surrounding anatomy. The honeycomb cells conform to body shape, and distribute weight over a wide area. Additionally, Stimulite Honeycomb stabilizes any position due to the axial buttressing of the honeycomb cell walls perpendicular to the anatomy.

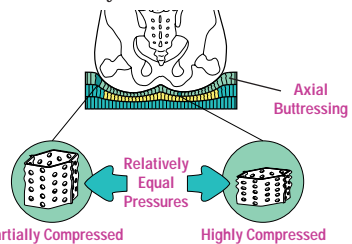

Uniform and Stable Load Distribution

Shear Force Reduction

The horizontal stiffness of a material determines both its shear resistance and stability. Too much stiffness increases shear forces, while too little decreases stability. Stimulite Honeycomb minimizes shearing without sacrificing stability. The honeycomb cells align themselves axially to the body and flex with movement to reduce shear forces against the skin.

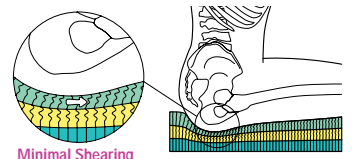

Diffuses Lateral Forces

Stimulite Silver Cushion



This affordable, comfortable cushion can be used in a wheelchair, on a recliner, in a car, plane or on a stadium seat and provides excellent stability. Made from a soft, flexible form of aerospace honeycomb, the lightweight Stimulite Silver Cushion distributes weight evenly to relieve pressure, conforms to body shape and flexes with movement to prevent shearing. At 1-1/2-inches thick and weighing only 1 pound, it's a cushion you can take with you wherever you go.

Like all Stimulite cushions, the Silver Cushion is naturally antibacterial, antifungal and odor resistant. *Machine washable* and *dryer safe*.

Stimulite Silver Cushion Sizes (width x depth)
16" x 16", 16" x 18", 18" x 16", 18" x 18"

To order, call 1.800.SUPRACOR or 1.800.787.7226.



800.787.7226
www.supracor.com